




## **BSE - CA Series Bedding Disposal Instrument**

## Bedding Collection Instrument



The BSE-CA series bedding collection table is specially used for the animal lab. It is used to prevent the staff from inhaling dust or polluted gas when dumping the cage bedding, and protect the safety of the staff.

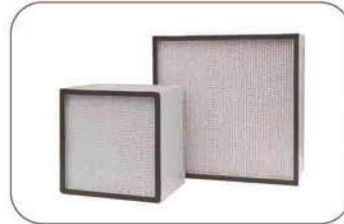

- Operating area for negative pressure, effectively adsorption of dust and odor, protect the safety of personnel.
- Primary effect/high efficiency double filtration, the polluted gas can be directly discharged into the room or connected with pipes to discharge the outdoor air after filtering by HEPA (99.995% high efficiency filter) to improve indoor air quality.
- SUS304 stainless steel dumping table, electrostatic spraying outer cover, easy to clean.
- The dumping port is horn shaped and inclined at a large angle, which is convenient for material collection.
- The dumping port is equipped with collision bar, cage box knock against collision bar, easy to remove silt pad material.
- There are large glass Windows on both sides of the equipment and lights on the top to increase light transmittance.
- Handrails are installed on both sides with casters at the bottom for easy movement.
- The bottom of the equipment is equipped with a large number of movable dirt barrels, which is convenient to collect waste mat bedding.



name	operator	Model	Overall dimensions ( L*W*H )	workspace dimensions ( L*W*H )	The collecting barrel
Bedding collection table	single	BSE-CA-A1000	1100 x 880 x 1950	1000 x 600	2

### Contact Information : Asian Biotechnology Inc

Unit # 001, Building A-1, Dubai Silicon Oasis DDP, Dubai PO Box No.342001 United Arab Emirates.

Email ID : [info@abidubai.com](mailto:info@abidubai.com) | Web : [www.abidubai.com](http://www.abidubai.com)

### Manufactured By: SHINVA MEDICAL INSTRUMENT CO., LTD

2009 Xinhua Road Zhoucun district Zibo City

Email ID : [zhouxin716@shinva.com](mailto:zhouxin716@shinva.com) | Web : [www.shinva.com](http://www.shinva.com)